



SuperPole™ for Angled Ceilings

An excellent option for rooms with angled ceilings. Works just like the award winning SuperPole™!

Features and Benefits

• Accommodates Various Ceiling Angles And Heights.

Angled ceiling plate can be mounted to any ceiling with a 45 to 90 degree angle. Fits standard 8ft (93"-99") ceilings or add the Uni-Fit Extender™ (sold separately) to increase ceiling heights up to 10ft (100"-120").

• Trouble-free Installation


Angled Ceiling Plate mounts securely to structural studs. Simply connect the SuperPole™ to the Angled Ceiling Plate and turn the jackscrew at the bottom of the SuperPole™ to create rock solid tension and support.

• Convert Any Standard SuperPole™


Any standard SuperPole™ can be converted into a SuperPole™ For Angled Ceilings by removing the standard top plate and replacing it with an angled ceiling top plate.

• Add On Kits For Evolving Needs.

The SuperPole™ For Angled Ceilings serves as a mounting structure for a complete range of options: The SuperBar™, Safety-Island™ and SuperTray™.

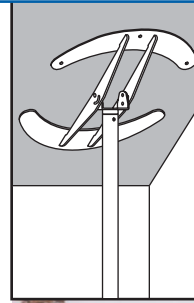

Jackscrew Installation

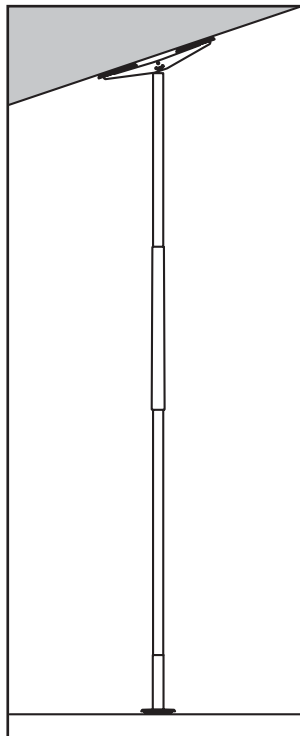

Technical Specifications:

SuperPole™	90" / 229 cm x 1.5" / 3.81 cm diameter. White powder coated steel.
Grip Length	20" / 51 cm. Washable, latex free grip.
Floor to Ceiling Range	93-99" / 236-253 cm (standard).
Expansion Mechanism	1-1/8" fine cut thread. Zinc plated steel jackscrew.
Top Plate	20" / 51 cm length. Top plate rubber will not mark ceiling.
Base Plate	5" / 13 cm diameter. Non-marking, base plate rubber ribbed for better gripping.
User Weight Capacity	300 lbs / 136 kg. Bariatric weight capacity 450 lbs / 205 kg.

For proper installation, top plate must be bolted into the angled ceiling. Not to be used with the SuperTrapeze™. Ensure ceiling structure and fastening means are capable of bearing intended support. For aggressive use, we suggest to secure the base plate to the floor.

SuperPole™ Accessories (each sold separately):

Uni-Fit Extender™	Increase floor to ceiling range from 100-120" / 254-305 cm. Not to be used with the SuperTrapeze™, Bariatric, & Ultra Models.
SuperPole™ Ultra	Longer pole accommodates floor to ceiling range from 121-140" / 307-356 cm. Not to be used with the SuperBar™, Safety Island™ or SuperTrapeze™.