



# SuperPole™ with SuperTrapeze™

An award winning trapeze with a unique offset handle that mounts to the SuperPole™. Provides superior assistance to achieve a seated position in bed.

## Features and Benefits

### ▪ Innovative.

Unlike conventional triangle trapeze handles, the SuperTrapeze™ handle simulates two rungs on a ladder. Allows for improved use of arm strength to “climb up” to a seated position.

### ▪ Optimizes Space.


Conventional trapeze floor stands can force the bed away from wall, reducing valuable floor space. SuperPole™ with SuperTrapeze™ installs beside the bed and requires less than 5" / 13 cm of floor space.

### ▪ Adjustable.

The trapeze support arm and nylon strap are height adjustable. For easy storage, simply hook the trapeze handle onto the SuperPole™.

MEDICARE BILLING CODE: E0940 (Freestanding trapeze).

Sold as SuperPole™ with SuperTrapeze™, or add SuperTrapeze™ to an existing SuperPole™.

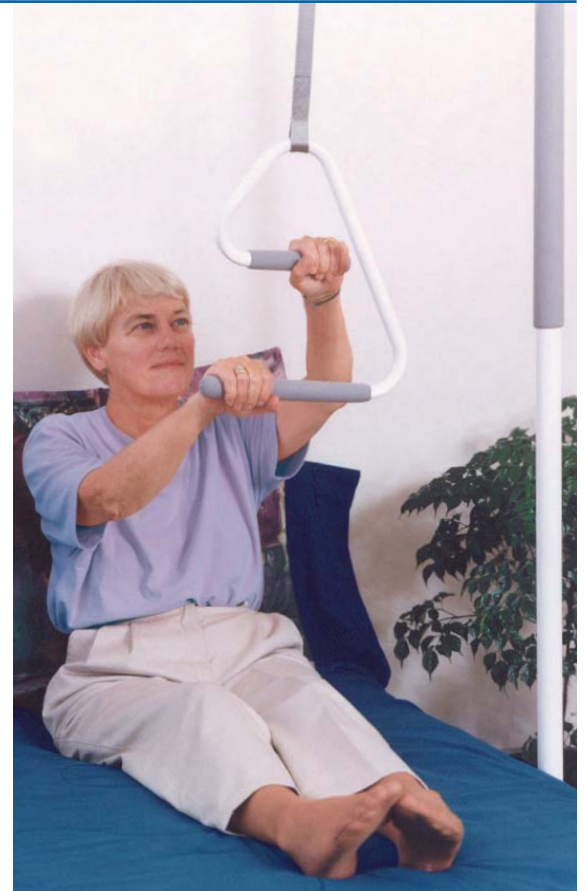

## Technical Specifications:



SuperTrapeze™	0.875" / 2.22 cm diameter rail. Strap range 42-70" / 107-178 cm from ceiling.
Trapeze Handle Grips	6.25" / 16 cm and 10.25" / 26 cm.
Support Arm	17" / 43 cm length. Recommended 84" / 213 cm from floor.
Strap Range	29" / 74 cm.
SuperPole™	90" / 229 cm x 1.5" / 3.81 cm diameter. White powder coated steel.
Grip Length	20" / 51 cm. Washable, latex free grip.
Floor to Ceiling Range	93-99" / 236-253 cm (standard).
Expansion Mechanism	1-1/8" fine cut thread. Zinc plated steel jackscrew.
Top Plate	20" / 51 cm length. Top plate rubber will not mark ceiling.
Base Plate	5" / 13 cm Ø. Non-marking, base plate rubber is ribbed for better gripping.
Weight Capacity	300 lbs / 136 kg.

No more than 150 lbs / 70 kg should be exerted on trapeze system (for repositioning and sitting up assistance only).

## SuperPole™ Accessories (each sold separately):

Uni-Fit Extender™	<b>Not to be used with the SuperTrapeze™, Bariatric, &amp; Ultra Models.</b>
SuperPole™ Ultra	<b>Not to be used with the SuperBar™, Safety Island™ or SuperTrapeze™.</b>
Ceiling Plate Extender™	Allows standard ceiling plate 20" / 51 cm length to span 30" / 76 cm with holes for ceiling attachment to 24" / 61 cm ceiling beam centers.