



QUICK MOUNT SUPERPOLE™ with TRI-FRAME CEILING MOUNT™

A permanent Tri-Frame mounted to the structural ceiling and concealed above acoustic tiled drop ceilings. Features a small docking station which extends below the drop ceiling and allows for a Quick Mount SuperPole™ to connect into place for transfer support.

TRI-FRAME CEILING MOUNT™

• Trouble-free Installation

Ceiling mast and threaded support rods mount into the structural ceiling and can be angled to avoid ducts, wiring and other obstructions. Offers a space efficient and economic alternative to lowering a ceiling to accommodate a SuperPole™

• Fits Virtually Any Drop Ceiling Range: 8"- 48"

• Docking Station

An inconspicuous docking station extends below the drop ceiling to allow for quick installation of a Quick Mount SuperPole™.

QUICK MOUNT SUPERPOLE™

• Quick Installation And Removal



Simply connect the Quick Mount SuperPole™ to the Tri-Frame dock and turn the jackscrew at the bottom of the pole to create rock solid tension. Supports bariatric patients up to 450lbs.

Convert Any Standard SuperPole

• Any standard SuperPole™ can be converted into a Quick Mount SuperPole™ by removing the standard SuperPole™ ceiling plate and replacing it with a quick-mount socket.

• Add On Kits For Evolving Needs.

The Quick Mount SuperPole™ serves as a mounting structure for a complete range of support options to reduce caregiver strain and improve patient mobility: The SuperBar™, Safety-Island™ and SuperTray™.

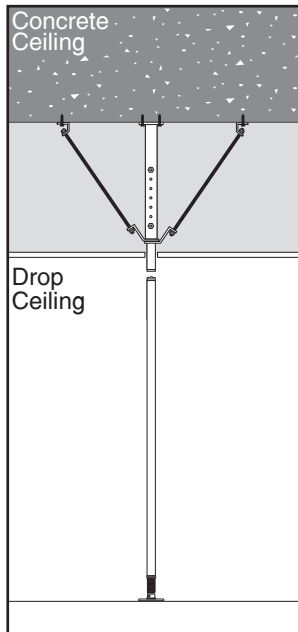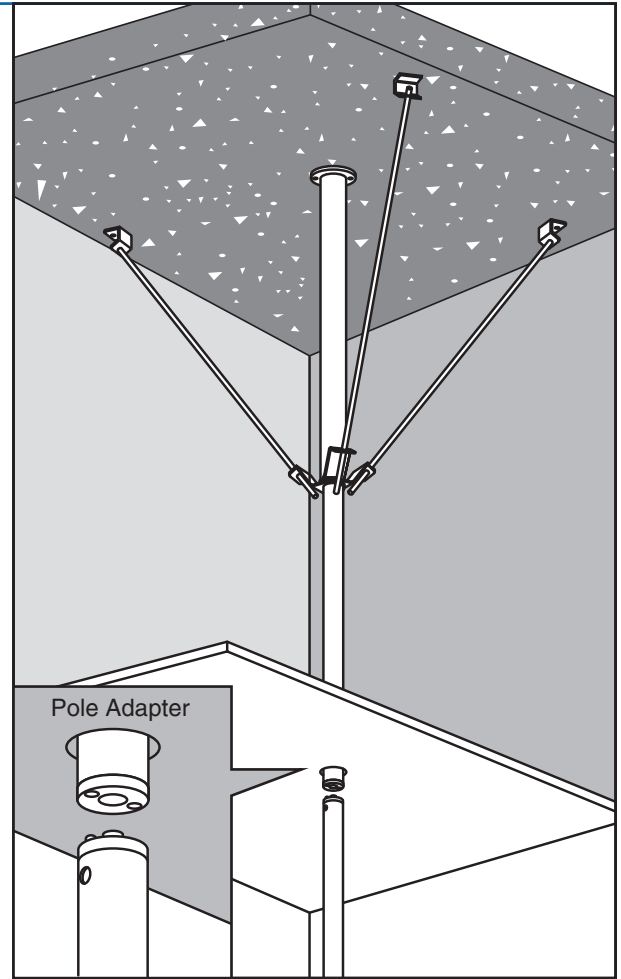

Technical Specifications:

SuperPole™	90" / 229 cm x 1.5" / 3.81 cm diameter. White powder coated steel.
Grip Length	20" / 51 cm. Washable, latex free grip.
Floor to Ceiling Range	93-99" / 236-253 cm (standard).
Expansion Mechanism	1-1/8" fine cut thread. Zinc plated steel jackscrew.
Base Plate	5" / 13 cm diameter. Non-marking, base plate rubber ribbed for better gripping.
User Weight Capacity	Bariatric weight capacity 450 lbs / 205 kg.