



e2™ Ceiling Mount Trapeze

A simple and economical alternative to floor base trapeze units. Features the innovative e2™ trapeze handle for increased leverage. Conveniently slides to the side of the bed when not in use.

Features and Benefits

• Innovative e2™ Trapeze Handle

Unlike conventional triangle trapeze handles, the e2™ trapeze handle simulates two rungs on a ladder allowing for better leverage when climbing up to a seated position.

• Slide-Away Convenience



Will not obstruct the view. When not in use, the e2™ trapeze handle conveniently slides to the edge of the bed.

• Space Efficient

Unlike bed mounted and floor based trapeze units that often require the bed to be away from the wall, The e2™ Ceiling Mount Trapeze maximizes valuable floor space.

• Installs Like A Grab Bar.

Mounts directly to structural studs of virtually any ceiling.

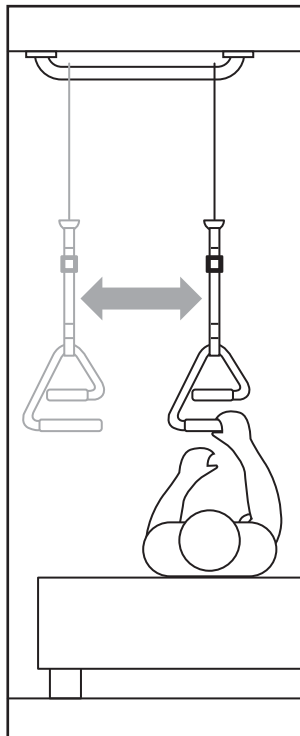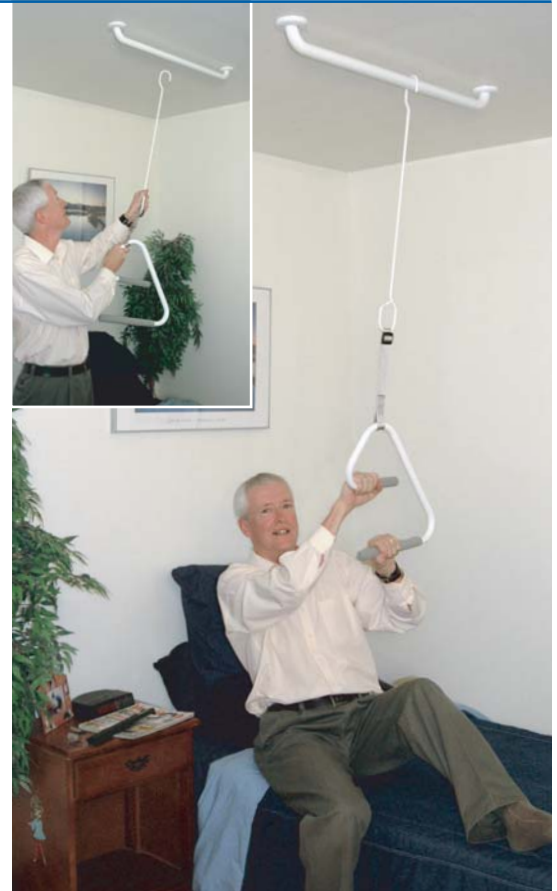

Technical Specifications:

e2™ Trapeze Handle	0.875" / 2.22 cm diameter rail. Grip lengths 6.25" / 16 cm and 10.25" / 26 cm.
Strap Range	3-29" / 7-74cm.
Ceiling Bar Diameter	1.25" / 3.2 cm.
Ceiling Bar Length	32" / 81 cm.
User Weight Capacity	300 lbs / 136 kg.

No more than 150 lbs / 70 kg should be exerted on trapeze system (for repositioning and sitting up assistance only). Ensure ceiling structure and fastening means are capable of bearing intended support.